Pets-Creating a bond

through the IoT-based digital healthcare service products

66

What is essential is invisible to the eye!

The Little Prince

- Antoine de Saint-Exupéry



















Little Cat

Our pets never TALK to us, even though they are very sick or hurt. So, For the never-verbal family member, 'the Little Cat' provide customized health-care solutions by gathering and analyzing data about pets through 'the Little Cat' IoT devices.

The Little Cat's all IoT products take a inspiration from The Little Prince story book. So, our all products are related to the characters and elements of The Little Prince.

the Little Prince



"We create a bond with pets

based on healthcare data"
















Medisand / Medipad

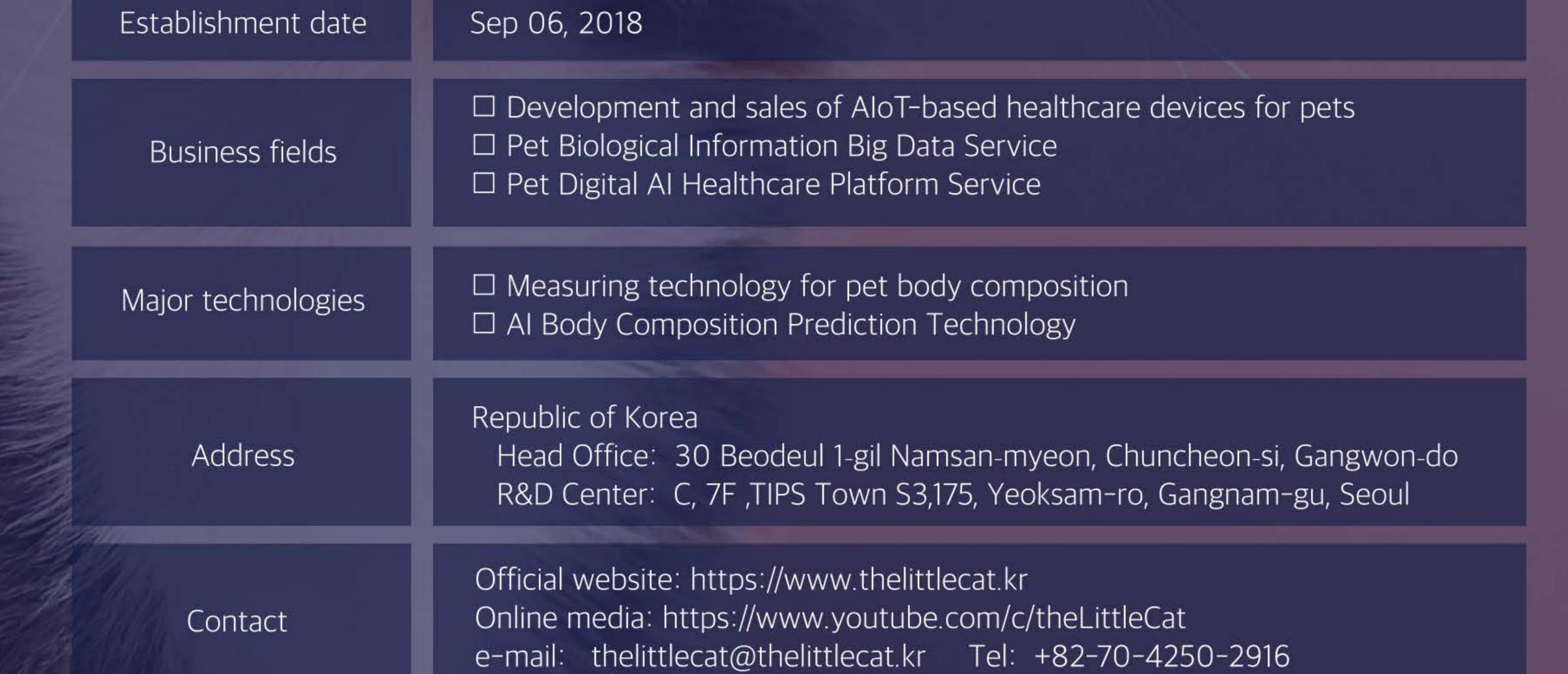




BIA(Bioelectrical Impedance Analysis) based on Cat Body Composition Monitor


Body composition monitor using cat's instinct


Taking advantage of the cat's instinct to enter the circle, it automatically measures body composition quickly and accurately through the paw pads when the paw enters the circle.

Monitor 8 body compositions with a dedicated app.

(Body water, Visceral fat, Subcutaneous fat, Protein, Muscle, Bone density, BMI, Weight)

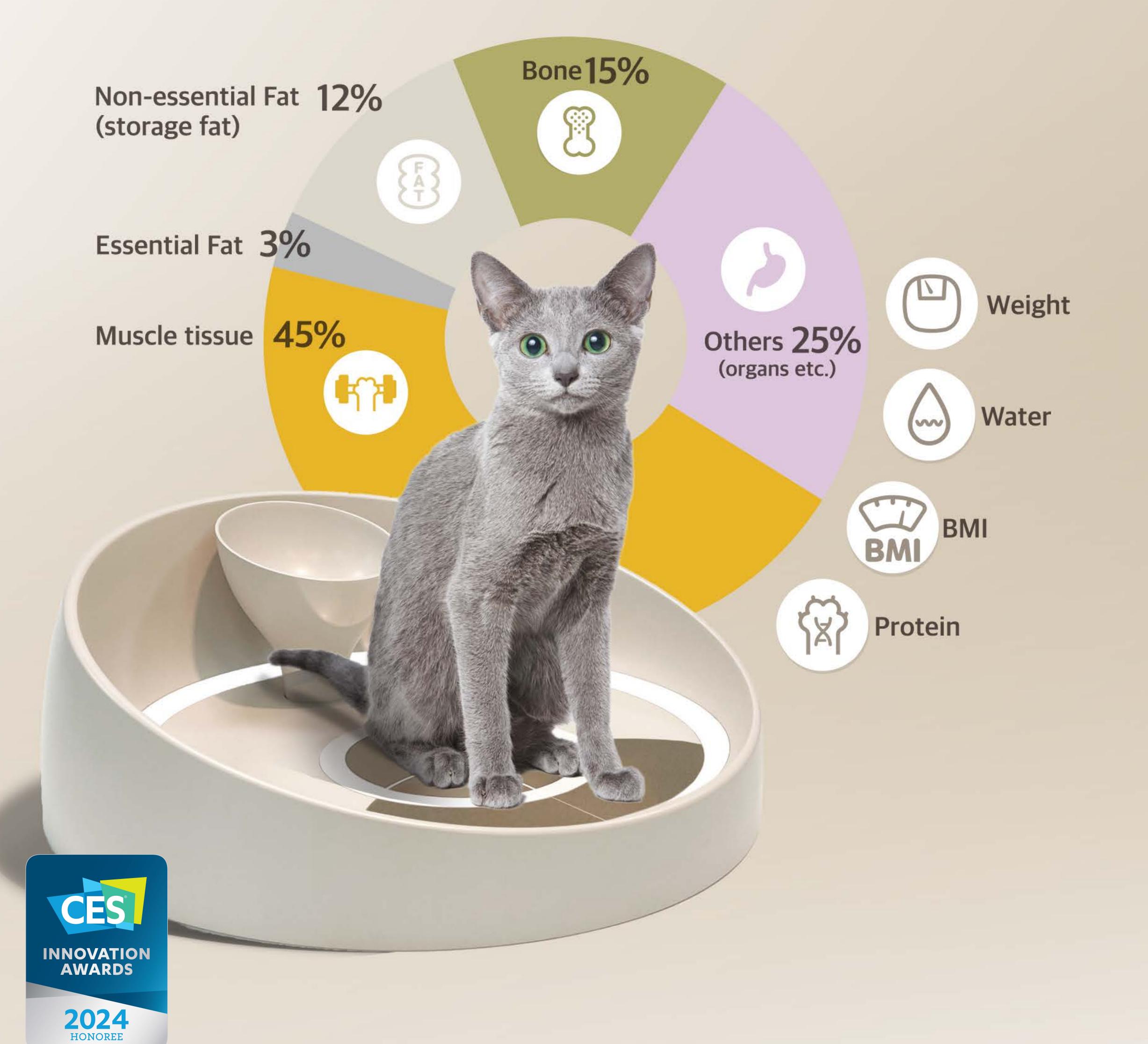

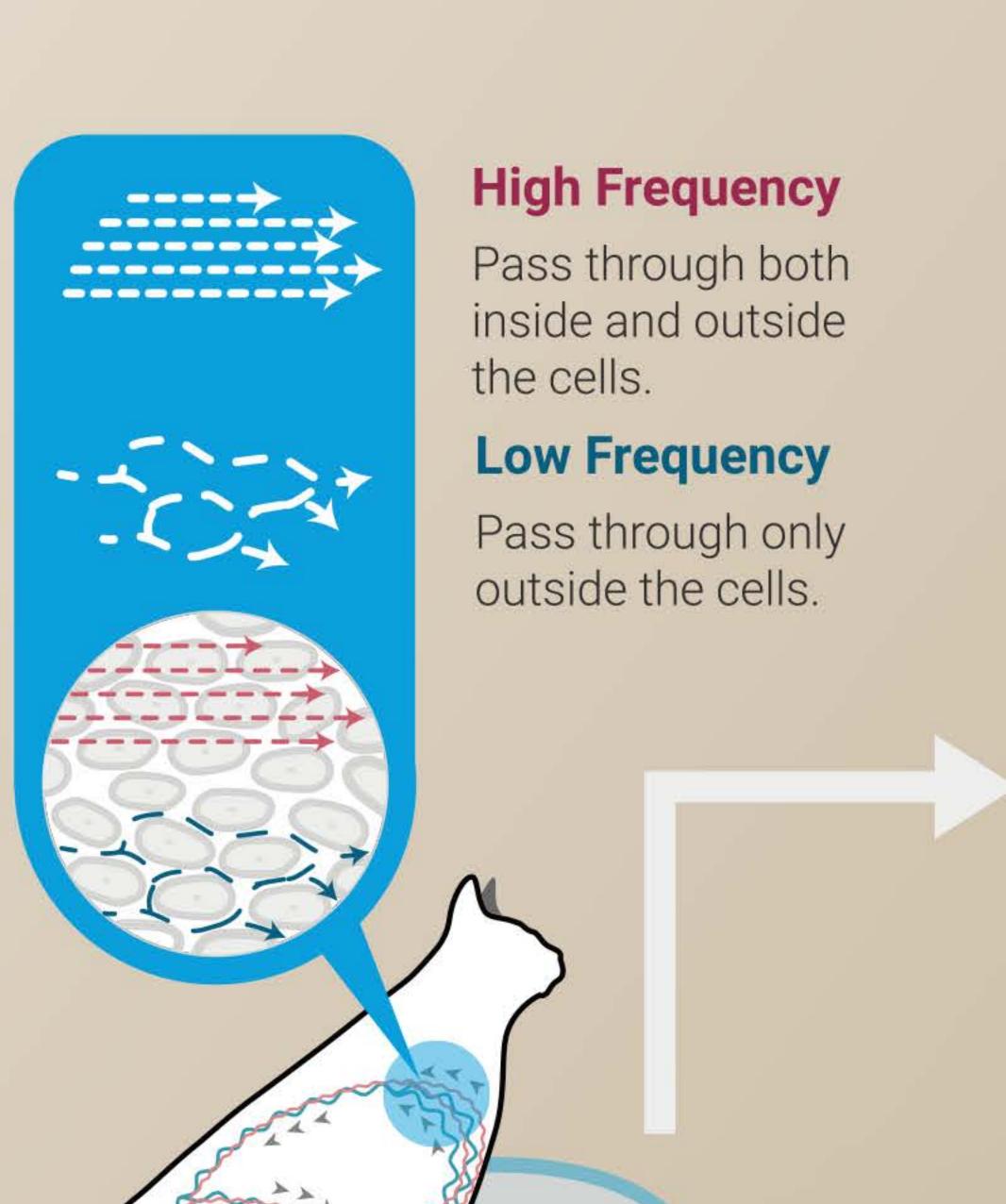

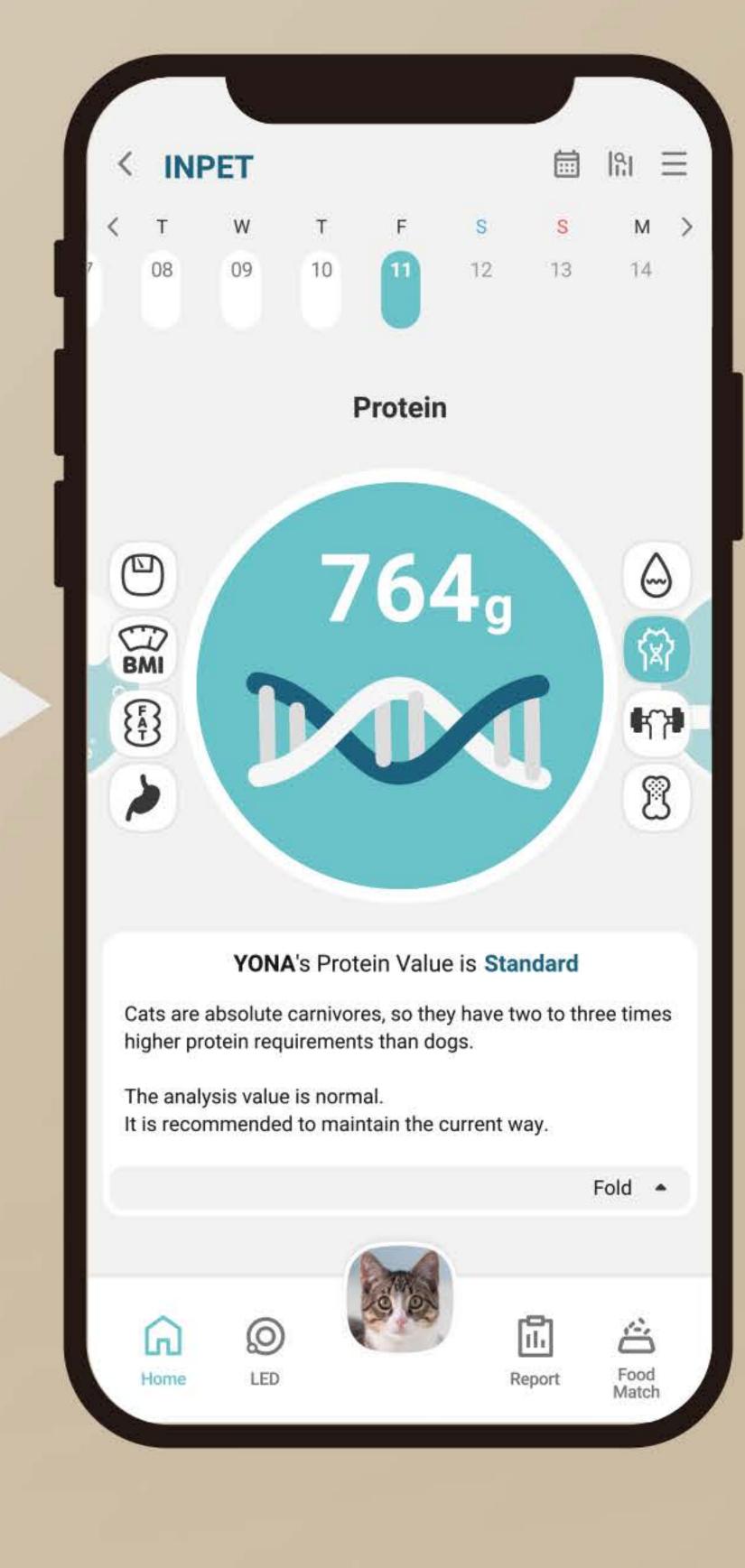



Weight BMI Fat VisceralFat

3.3

40%

70%

Water Protein Muscle Bone



Feed recommendations based on body composition



Vital signs and telemedicine alarms

Reduce healthcare costs by predicting signs of illness through changes in body composition, supporting telemedicine consultations, and preventing disease



Premium discounts

Reduced premiums with regular body composition reports and real-time monitoring



Cat Smart Wheel

IoT-based indoor exercise machine linked with application Designed to stimulate the dynamic visual acuity of cats



Easy to useYou can easily check the cat exercise data with the app.



Bluetooth connection

Bluetooth Connectivity to Measure and Adjust Data Value



LED light controlYou can adjust the LED light color and light position with the app.



Exercise measurement

You can check the calorie consumption according to the amount of exercise.




Data accumulation
All your cat's exercise data is stored daily.



Cat information
Name, Birthday, BCS(Body Condition Score)





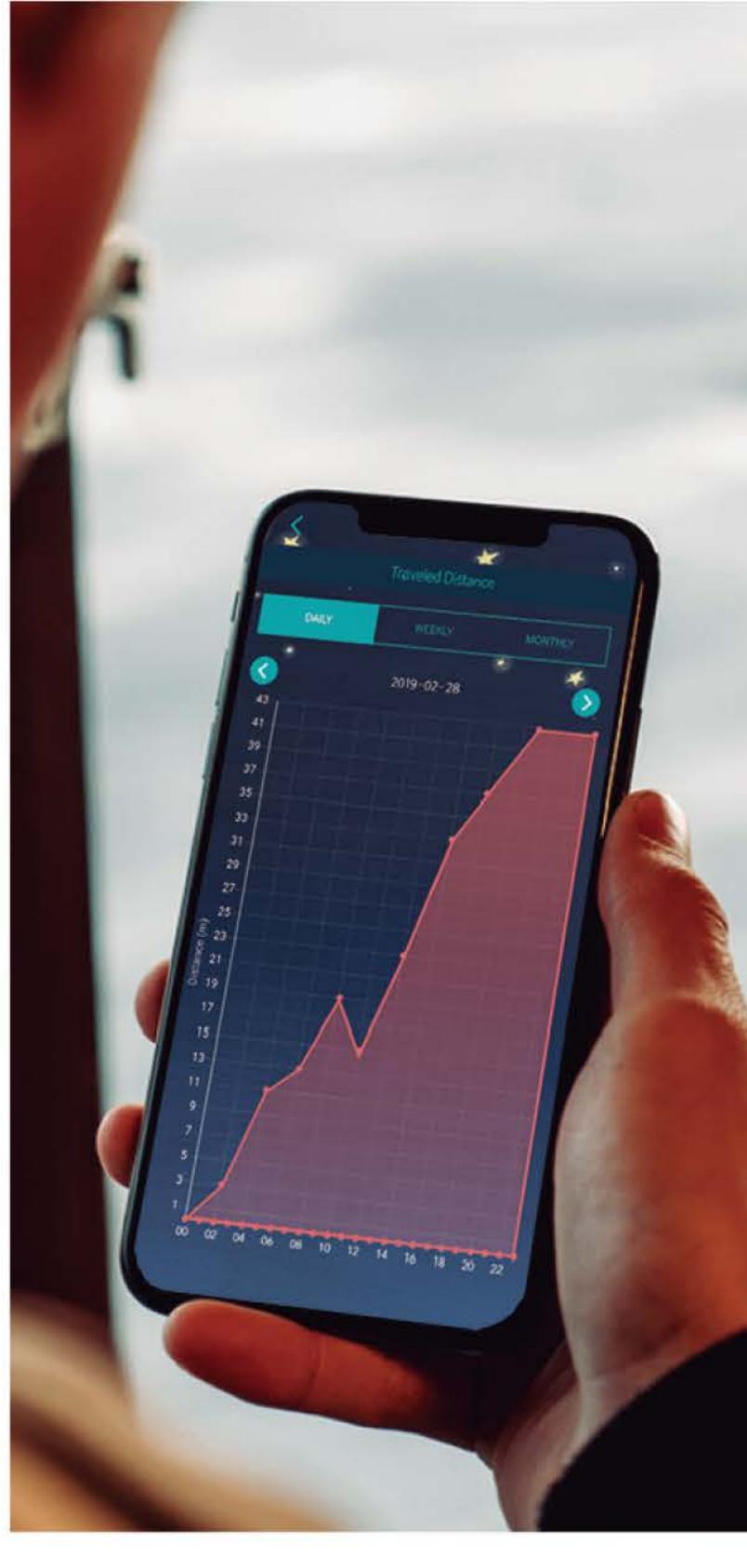





Dog Health Checking Pee Pads

A pee pad that checks the dog's kidney health in daily life



Pad pattern color changes by kidney disease

Pink - Alkali reaction Blue - Blood reaction



Absolutely chemical-free

No chemical coagulants No chemical absorbents No chemical odorants




Check you dog's kidney health with our app

App: Medipad You can check your dog's kidney health with MEDIPAD app.




PE waterproof film

Excellent waterproofing PE film, thez floor does not get wet.

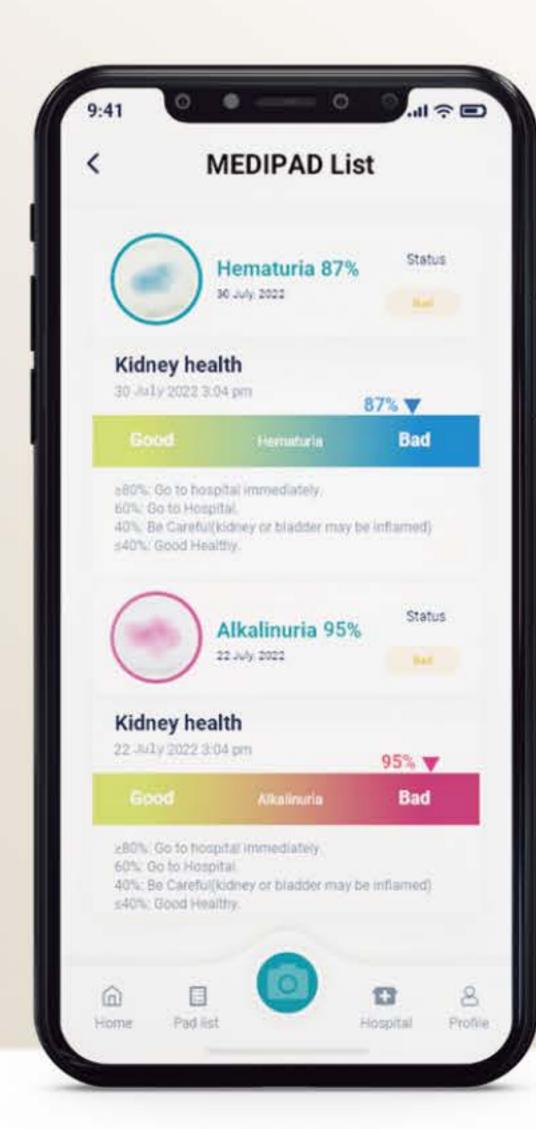




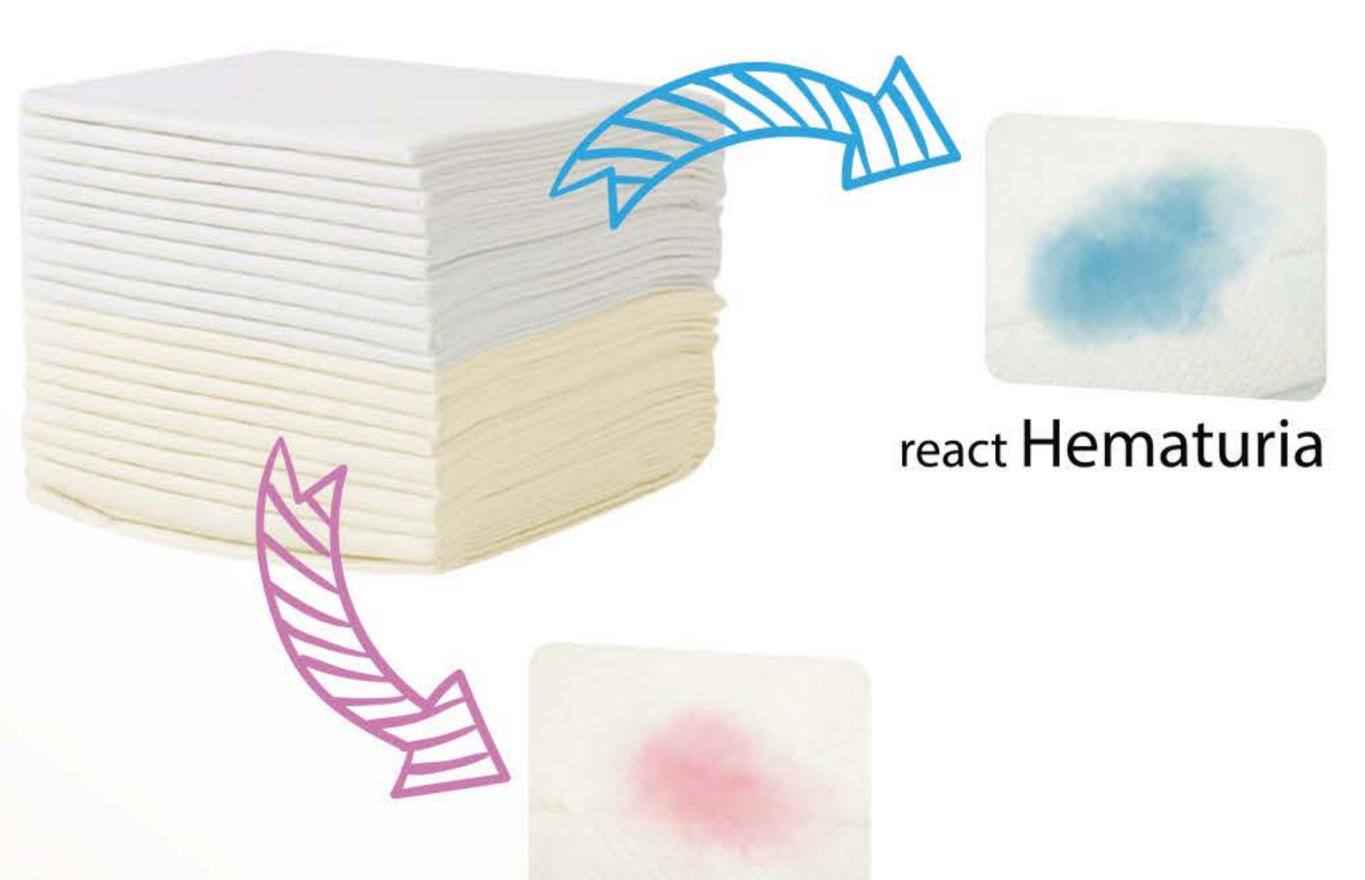





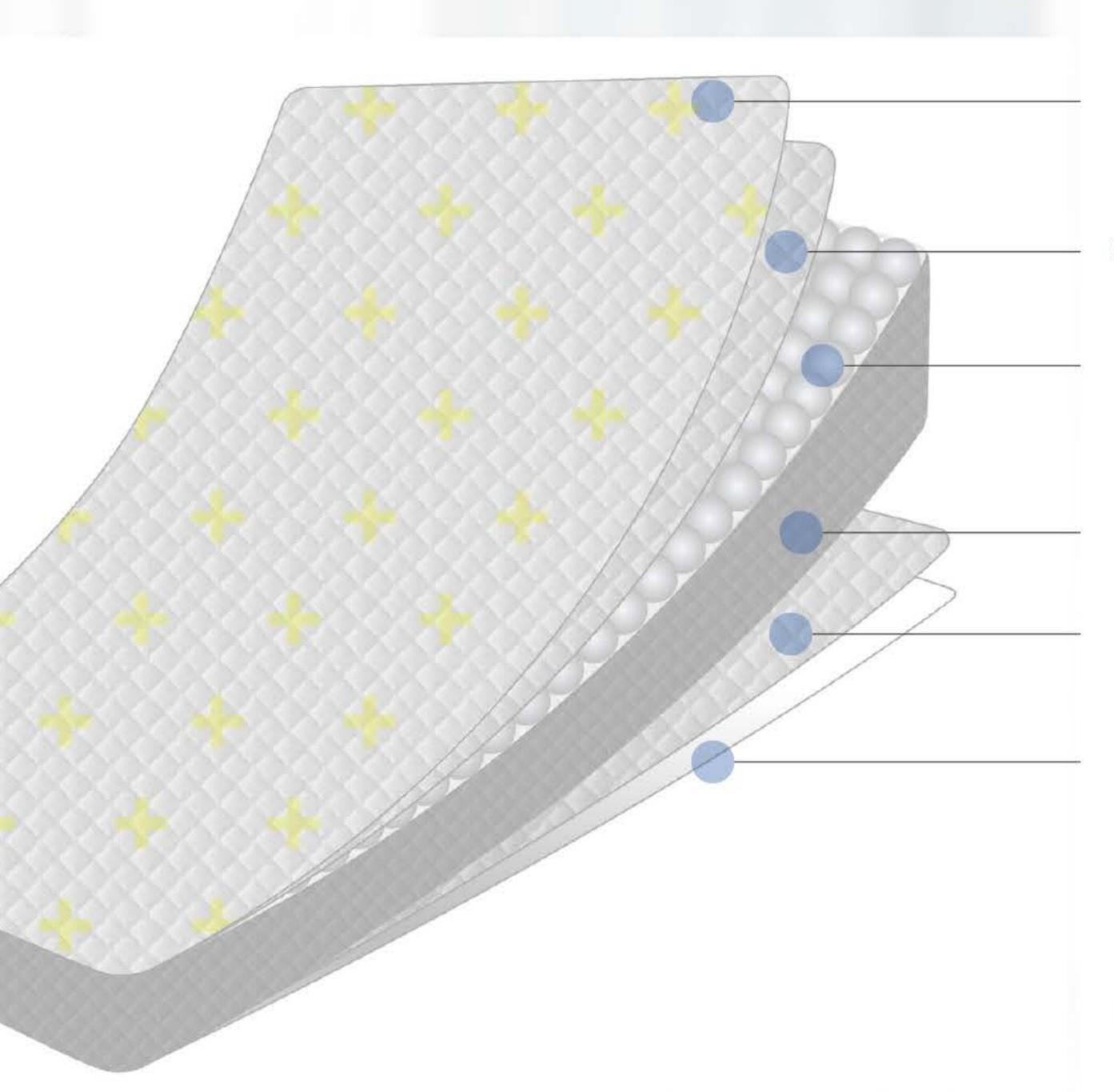

Alkaline Urine Reaction Skin Cover

Square embossing pattern prevents urine from spreading to both sides

Compressed Square Embossed Paper

Super Absorbent Polymer(SAP)
It contains 3 to 4 times more than other products, so it absorbs quickly and has a lot of

100% Pure Wood Pulp

Compressed Square Embossed Paper

Waterproof Film. The floor does not get wet.



Model	MEDIPAD(MEDIPAD-A4050, MEDIPAD-A6076)
Size	• MEDIPAD A4050: 400mm x 500mm • MEDIPAD A6076: 600mm x 760mm
Material	• 6 Layers: PP + Paper + SAP + Pulp + Paper + PE • Package: PE Film
Package	 MEDIPAD A4050 100Pcs (Normal 90 + Test Pad 10) /Pack, Pack: 280mm x 140mm x 320mm, Carton Box:440mm x 320mm x 320mm(3 Packs) MEDIPAD A6076 25Pcs (Normal 22 + Test Pad 3)/Pack, Pack: 205mm x 160mm x 250mm, Carton Box:420mm x 260mm x 360mm(4 Packs)
Feature	 React Alkalinuria Produced as an OEM in China. MediPad is a simple urine test for dogs that checks the degree of alkalization of their urine. Medi pad is not a kit for diagnosing diseases, nor is it for medical use. MediPad was designed to check that dogs' urine is alkaline in the early stages due to infection with bacteria or viruses of lower urinary tract diseases (urinary tract infections, metabolic disorders, etc.)
Notice	• HS CODE: 9619.00




Depending on the degree of alkali in the dog's urine, the cross-shaped yellow color changes to a pink color.

Cat Health Checking Sand

Sand to check your cat's kidney health in everyday life



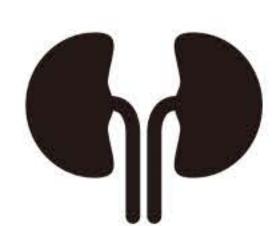

Bead color changes by kidney disease

Pink - Alkali reaction Blue - Blood reaction




Absolutely chemical-free

No chemical coagulants No chemical absorbents No chemical odorants



Check you cat's kidney health with our app

App: Medisand

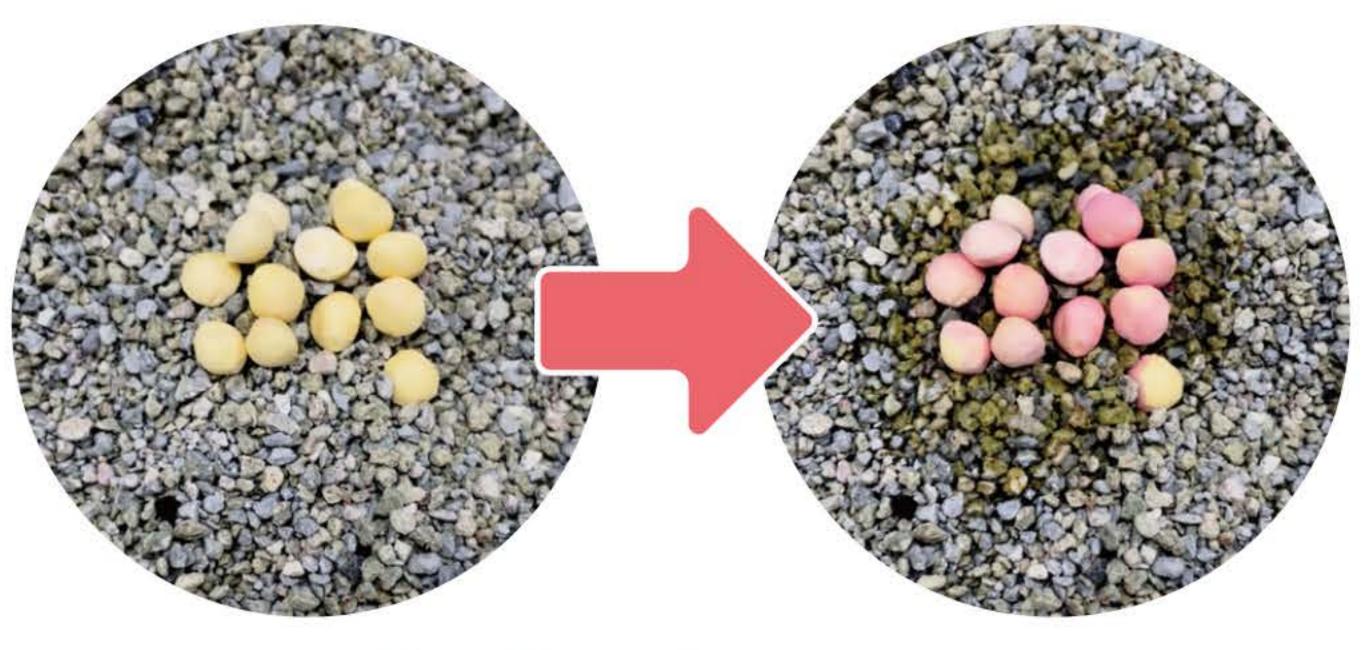

Alkali-pink reaction



Blood - blue reaction

Record Data

kidney health.


Description

Mobile application-linked


enables you to check their

Depending on a cat's kidney health, the diagnostic bead in the sand respond to color, which enables the owners to check the level of their cat's health more accurately via a mobile application.



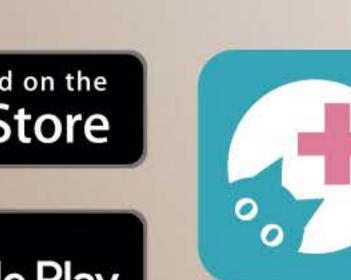




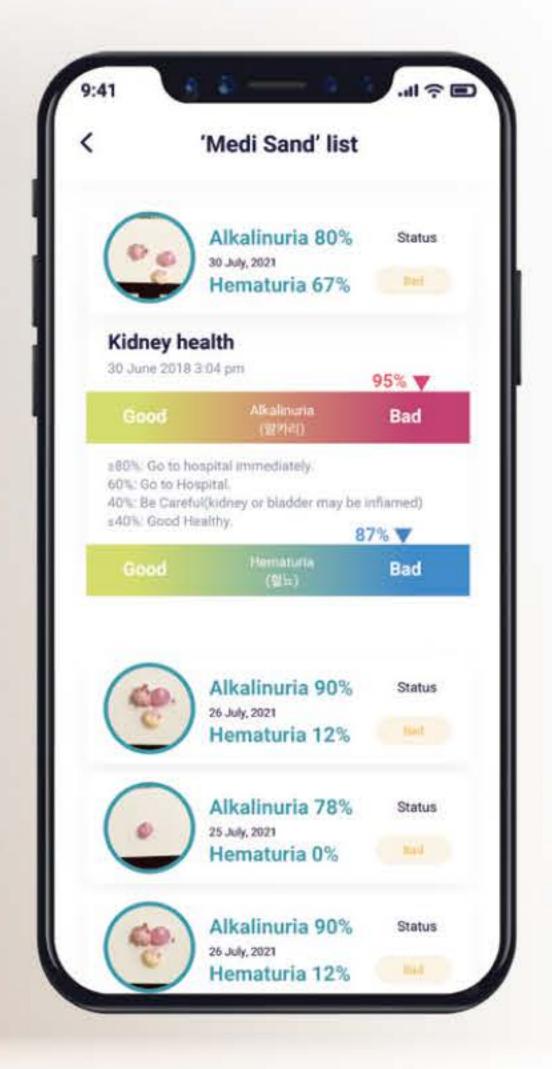



Model	MEDISAND(MEDISAND-A200, MEDISAND-H200)
Size	 Product Size: 125mm x 32mm x 190mm(Inner Box Size), 200g/ea Package(Carton Box) Size: 370mm x 270mm x 400mm, 40Ea
Material	Bead: Potato & Corn starch Package: Zipper pack(PE+ Aluminium)
Feature	 Color changes in hematuria and alkaline urine MEDISAND is a simple urine test for cats that checks the degree of alkalization and hematuria of cats' urine, not medical tests. The final diagnosis should be made by a veterinarian's opinion and examination at the hospital.
Notice	• HS CODE: 3824.99

How to Use

- Clean the cat bathroom you are using and put new litter on it.
- 2. Sprinkle MediSand evenly on the cat litter.
- Check if there is a change in the color of MediSand that came into contact with the urine after the cat urinated.
- If there is a color change of MediSand, observe it repeatedly for at least one day.

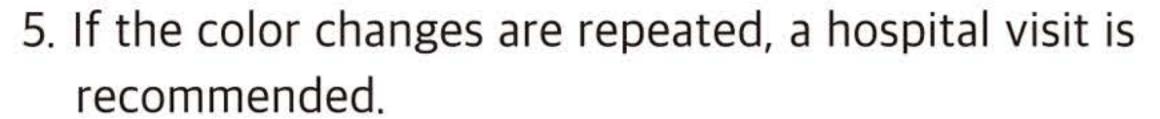

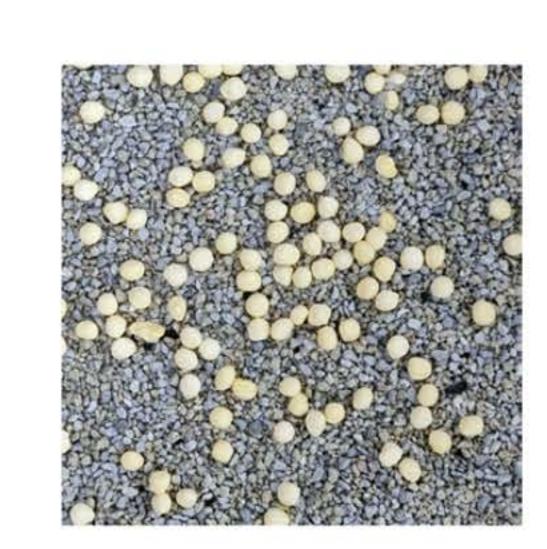


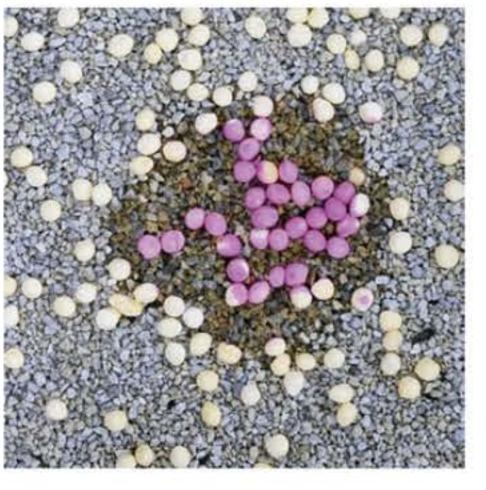




All grown-ups were once children


-The Little Prince

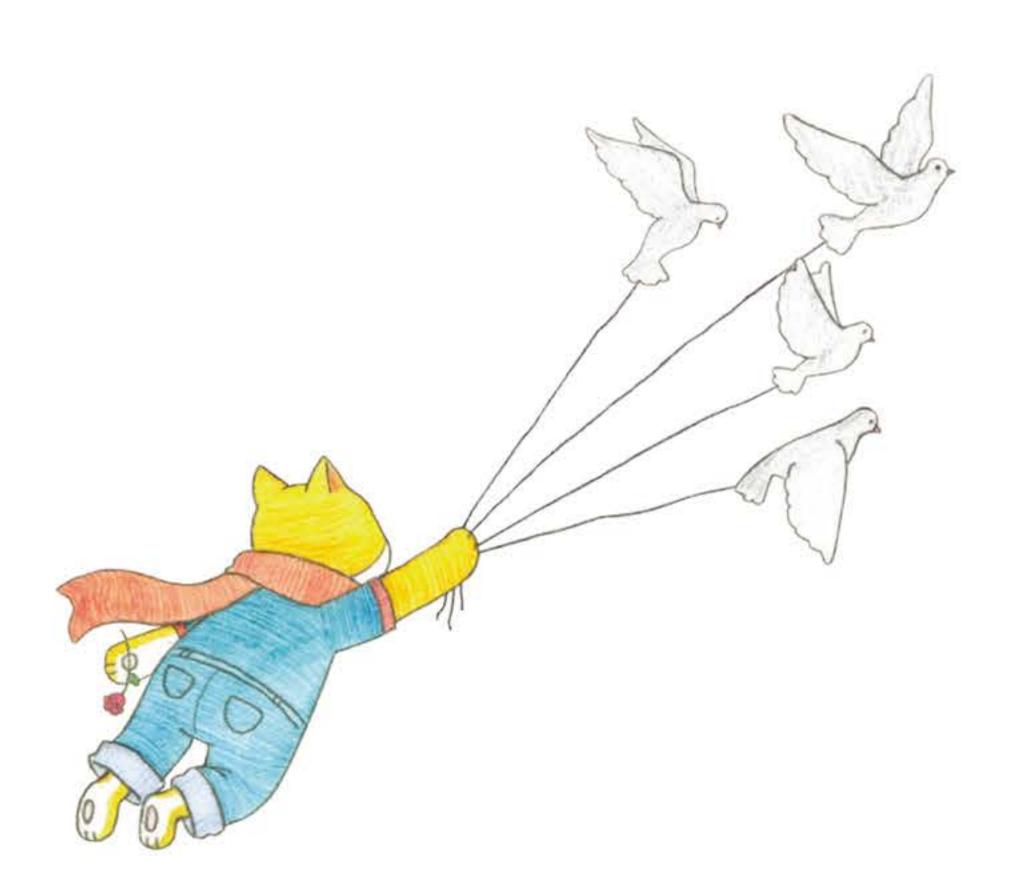



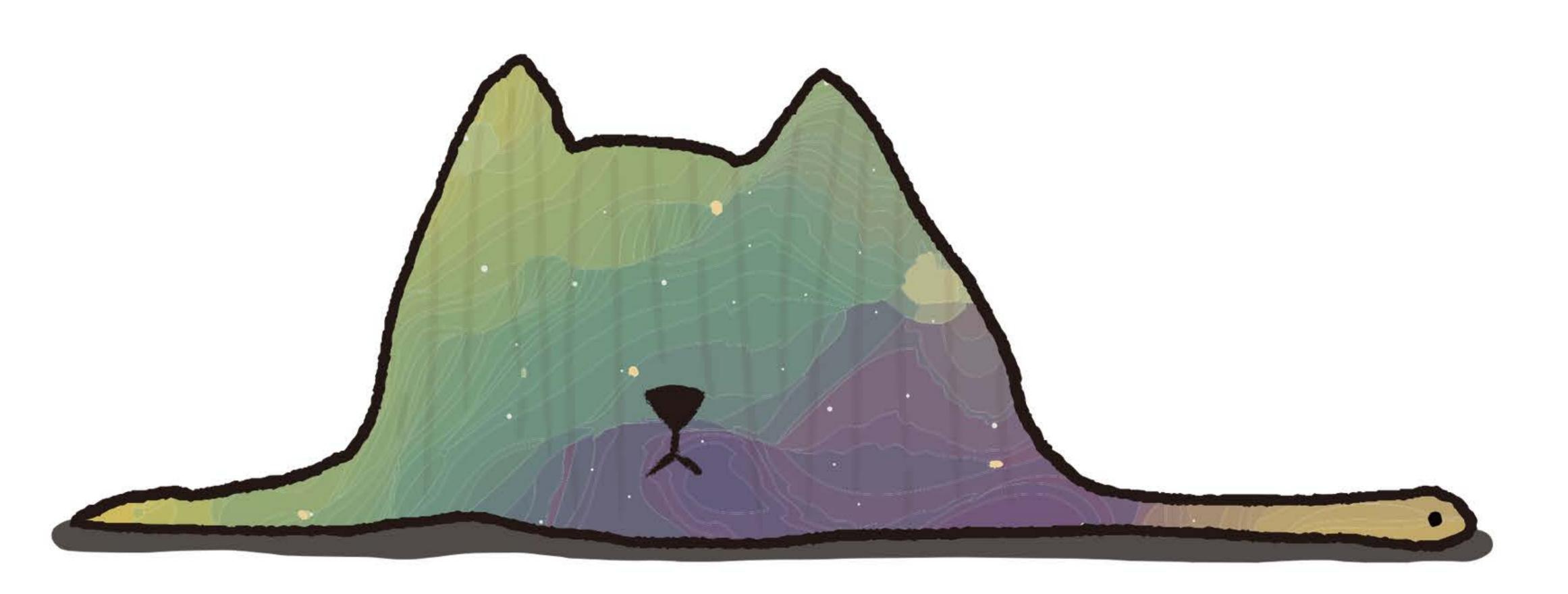

It was a picture of a boa constrictor digesting a ?



"Anything essential is invisible to the eyes."

-The Little Prince

MADE IN KOREA

Certified by the Korean Government,
National brand SME co-brand [Brand K]
Certified Brand [the Little Cat]



Pets health care service







the Little Cat

HO: Douzone BI center, 130 Beodeul 1-gil,Namsan-myeon, Chuncheon-si, Gangwon-do Rep. of KOREA Lab: B05,2F,2dong,Startup Campus,20, Pangyo-ro 289beon-gil, Bundang-gu, Seongnam-si, Gyeonggi-do

Tel +82-70-4250-2916 Fax +82-504-288-2916

the Little Cat. Co.,Ltd. www.thelittlecat.kr ©.2023.09